Health Kit

INSTRUCTION MANUAL

Medicine Spoon & Medicine Dropper

- Child should sit up while medicine is being dispensed.
- Dispense slowly.
- Thoroughly clean before and after each use.

CAUTION

Keep Medications out of your child's reach. It is recommended to use a child resistant top on all medicine bottles and keep them stored in an area that is not accessible by a child.

Medicine Spoon

To Use:

While child is sitting upright, hold the spoon upside down. Pour the medicine into the handle. Use the graduation to measure the proper amount. Level the spoon out and place into child's mouth. Allow child to swallow naturally. Do not pour the medicine into child's mouth, as this may result in choking.

To Wash:

Wash thoroughly in warm, soapy water before and after each use. The Medicine spoon is also dishwasher safe.

Medicine Dropper

To Use:

Squeeze the bulb at the end of the dropper and insert the tip into the medicine bottle. Releasing the bulb fills the tube with medicine. Fill bulb until the desired dosage marking. To eliminate excess, squeeze the bulb. While child is in an upright position, administer the medicine by inserting the tip of the dropper into the child's mouth and gently squeeze the bulb.

To Clean:

Before and after each use, place the medicine dropper in warm soapy water. Squeeze and release the bulb to flush the interior. To rinse, repeat the process in clear water. The Medicine Dropper is also dishwasher safe.

Nasal Aspirator

To Use:

While child is sitting upright, squeeze bulb to empty the air out. Gently insert the tip of the Aspirator into child's nostril. Slowly release the bulb to clear out child's passages. Repeat in other nostril.

To Clean:

Before and after each use, place in warm soapy water. Squeeze and release the bulb to flush interior. To rinse, repeat process in clear water. Be sure all water is drained during drying.

CAUTION

Never insert the aspirator too deep into the nostril. Always dry after each use. Not for use in child's ears or mouth. Keep away from heat and bright light.

Meets
CPSC Safety
Requirements

Thermometer

WARNINGS BEFORE USING DIGITAL THERMOMETER PLEASE READ CAREFULLY

- This is **NOT** an ear thermometer. **NEVER** put probe into ear.
- Take temperate measurements Underarm, Orally, or Rectally.
- **ALWAYS** supervise your children when using thermometer.
- **NEVER** allow your child to use it alone.
- **NEVER** allow children to play with this product. It is not a toy.

IMPORTANT INFORMATION

- Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.
- **DO NOT** allow child to walk or run during temperature taking. Probe could present a hazard.
- High, prolonged fever requires medical attention, especially for young children. Please contact your physician.
- Fever chart included in manual.


CARE & MAINTENANCE

- Clean the thermometer tip and probe using isopropyl alcohol or soap and water.
- Store the unit in the protective case while not in use.
- Do not store or keep in direct sunlight, or in high temperatures.
- Keep battery away from children
- To avoid damage to probe, do not bend the thermometer tip.
- Do not open the thermometer except to replace battery.
- Do not drop thermometer, as this may ruin the unit.

DIRECTIONS

RECTAL USE


- Commonly used for babies, young children, or when it is difficult to take an oral or underarm temperature.
- Lubricate probe tip with a water-soluble lubricant or petroleum jelly.
- Gently insert the probe tip no more than 1.3cm (1/2inch) into the rectum. Do not force the tip into the rectum if resistance is encountered.
- Clean the thermometer as recommended above.
- The peak temperature should be reached in approximately 15 seconds.
- Read the digital display and subtract 0.5°C (0.9°F) for the final temperature reading.

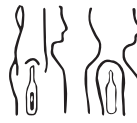UNDERARM (AXILLARY) USE

(For children 4 years and younger)

- Wipe underarm with a dry towel
- Place probe tip in armpit and bring arm down against body. The thermometer should be pointing upward in the same direction as your body if possible.
- The time required for an accurate underarm temperature may be 25 seconds or longer depending on placement of the probe and closure of the arm against the body.
- Read the digital display and add 0.5°C (0.9°F) for final temperature reading.

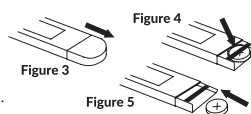CHANGING THE BATTERY

- Replace battery when " " appears in the lower right corner of LCD display.
- Pull battery cover off as shown in Figure 3.
- Gently pull out plastic circuit board with battery chamber approximately 1 cm (slightly less than 1/2") (See Figure 4)
- Use pointed object such as a pen to remove old battery. Discard battery lawfully. Replace with new 1.5V DC button type LR41 or SR41,UCC392, or equivalent. Be sure battery is installed with "+" polarity facing up. (See Figure 5)
- Slide battery chamber back into place and attach cover.

WARNING: KEEP BATTERY OUT OF CHILD'S REACH. SWALLOWING BATTERY COULD BE FATAL. BATTERY SHOULD NOT BE CHARGED OR PLACED INTO EXTREME HEAT AS IT MAY RUPTURE.

This is a sensitive electronic device and may be damaged if not used properly under parental supervision.

- Do not expose the thermometer to extreme temperatures, humidity, direct sunlight. Avoid sudden impact and dust.
- Keep battery out of child's reach. Swallowing a battery could be fatal.
- Batter should not be charged or placed into extreme heat, as it may explode.
- Used batteries should be disposed of properly.

PRECAUTION

The performance of the device may be degraded should one or more of the following occur:

- Operation outside the manufacture's stated temperature and humidity range. (41°F ~ 104°F)
- Storage outside the manufacturer's stated temperature and humidity range. (-4°F ~ 131°F)
- Mechanical shock (for example, drop test)
- Patient temperature is below ambient temperature.

Portable and mobile RF communications can affect the device. The device needs special precautions regarding EMC according to EMC information provided in the accompany documents.


WHAT IS A "NORMAL" TEMPERATURE?

Temperature readings vary in different individuals. A normal temperature refers to each individual's average temperature while they are well. Normal average temperatures range from 36.1°C to 37.7°C (97°F to 100°F) when taken orally. See "RECTAL USE" and "UNDERARM (AXILLARY) USE" for the oral equivalent of a rectal or underarm reading.

TEMPERATURE METHOD	NORMAL TEMPERATURE	FEVER RANGE
Oral	97.0°F ~ 99.6°F	99.6°F +
Rectal	97.3°F ~ 100.4°F	100.4°F +
Underarm	95.4°F ~ 98.6°F	98.6°F +


ORAL USE (For children 4 years and older)

- Do not drink hot or cold fluids, exercise, or perform other activities that will raise or lower temperature readings when compared to your normal, average temperature. The mouth should remain closed up to 5 minutes before attempting a reading.
- Place the probe tip well under the tongue, (see diagram) as indicated by the "✓" symbol. Instruct children not to bite down on the probe.
- The peak temperature should be reached in approximately 25 seconds. Opening the mouth could result in a longer time for a reading.
- When peak temperature is reached, the °C or °F symbol will stop flashing, and the thermometer will beep repeatedly. **Please note: It is not necessary to hear the beep to ensure the peak temperature has been reached. When the °C or °F symbol stops flashing and the display reads a consistent temperature, the peak temperature has been reached.**
- Read and record the temperature and time for reference. The reading will automatically be stored in the memory.
- Push the ON - OFF button to turn the unit off. If you do not turn the unit off, it will shut off automatically in approximately 10 minutes.

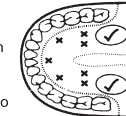TROUBLESHOOTING

Error	Problem	Solution
Lo	Temperature taken is lower than 32.0°C (90.0°F)	Turn off, wait one minute and take a new temperature via close contact and sufficient rest
Hi	Temperature taken is higher than 42.9°C (109.9°F)	
Err	The system is not functioning properly.	Unload the battery, wait for 1 minute and repower it. If the message reappears, contact the retailer for service.
	Dead battery: Battery icon is flashing, can't be measurable.	Replace the battery.

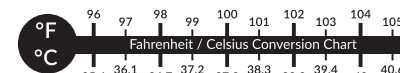Visit us at CribsforsKids.org

CAUTION

These products are not toys.
Do Not permit children to play with them.
When not in use, keep out of reach of children.

Nail Clippers

To Use:

1. Wash before initial use.
2. Insert the edge of baby's fingernail between the curved blades.
3. To clip, squeeze the clipper's handle and base between thumb and forefinger.

Note:

Do not clip too close to baby's skin.
We recommend clipping baby's nails while he or she is sleeping.

To Clean:

Sterilize with isopropyl rubbing alcohol and rinse in warm water. Wipe dry.

Toothbrush

To Use:

1. Begin use once all primary teeth are in or when recommended by your dentist.
2. Apply a pea-sized amount of toothpaste to the brush and have your child brush for 2-3 minutes using an up and down motion.
3. Continue to help and supervise your child during brushing - do not allow child to swallow toothpaste.

Note:

Toothbrushes should be replaced every 2-3 months.

To Clean:

Rinse each brush thoroughly before first use and after each use.

IMPORTANT

Questions about your child's oral health should be answered by your family dentist.

CAUTION

Adult supervision required.
Do not allow child to run or walk while carrying or using toothbrush.
Do not let children play with toothbrushes.

Korbex International LLC
3422 Old Capitol Trail, Suite 384
Wilmington, DE 19808 USA
1-888-608-7726
www.korbexbaby.com

Item Code: HK3001
Style# ETC
Lot# 052623
Date of Prod.

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